


EXIT-PRO[®] AG

SoftTouch

Antiseptic silver dressing for exit-site with cut and hole set on hypoallergenic adhesive

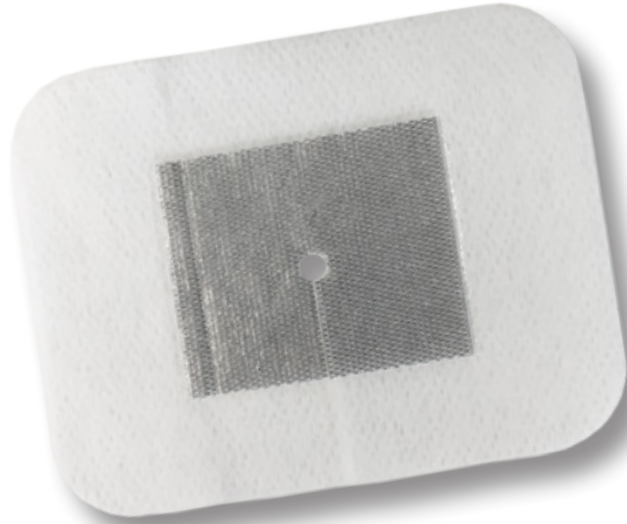

Description

Exit-site dressing with pre-cut silver ion-releasing anti-adherent film for the catheter exit, applied on a transpiring non-woven patch coated with silicone-based hypoallergenic adhesive gel with **SoftTouch** technology.

Intended use

Complete dressing of the insertion points of intravascular devices. It does not require the use of secondary dressings. Compatible with all types of catheters of standard use.

Indication

The dressing performs a double action:

- the silver ion-releasing film has an immediate and long-lasting antimicrobial action, helping to reduce the risk of local infections of exit-site and tunnel as well as the sepsis associated with the use of intravascular devices.
- the **SoftTouch** adhesive film preserves the skin from mechanical damage caused by continuous removal of the patch (as a consequence of repeated treatments), helping tissue repair.

Optimal solution for patients experiencing dry, moist or scaly skin widespread in patients undergoing, for example, dialysis treatment.

Available Models

CODE	MODEL	SILVER PAD	NW PLASTER	PCS/PACK	PACKS/BOX
S0810/AG	Cut and Hole ø 4 mm	5 x 4 cm	10 x 8 cm	1	50

⚠ Do not use silver dressings on patients sensitive to silver



Made in Italy by

EMODIAL

SOFTTOUCH LINE

