



Made in Italy by



EXIT-PRO[®] AG
SoftTouch



Catheter Dressing




www.valiantmedical.co.uk

+44(0)330 303 9090

info@valiantmedical.co.uk

VM-EXITPRO AGST-BROCHURE-04/21

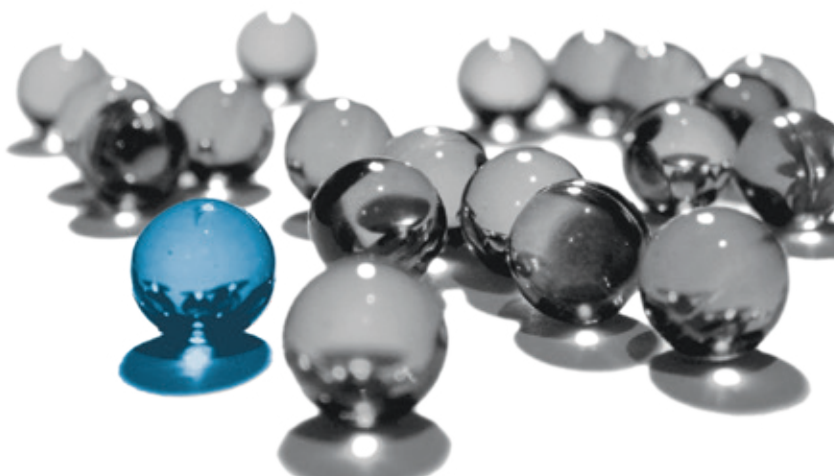

DESCRIPTION

EXIT-PRO AG SoftTouch is an exit-site dressing with pre-cut silver ion-releasing anti-adherent film for the catheter exit, applied on a transpiring non-woven patch coated with silicone-based hypoallergenic adhesive gel with "SoftTouch" technology.


INTENDED USE

Complete dressing of the insertion points of intravascular devices. It does not require the use of secondary dressings. Compatible with all types of catheters of standard use.

INDICATION


The dressing performs a double action:

- the silver ion-releasing film has an immediate and long-lasting antimicrobial action, helping to reduce the risk of local infections of exit-site and tunnel as well as the sepsis associated with the use of intravascular devices.
- the "SoftTouch" adhesive film preserves the skin from mechanical damage caused by continuous removal of the patch (as a consequence of repeated treatments), helping tissue repair. Optimal solution for patients experiencing dry, moist or scaly skin widespread in patients undergoing, for example, dialysis treatment and oncology treatment.

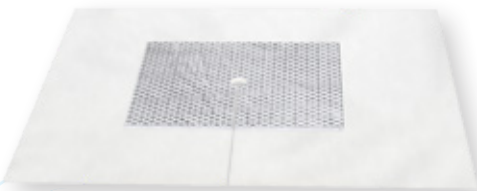

The softness of SoftTouch adhesive is combined with the antibacterial and preventive action of silver.

Ideal for chronic patients or for patients with delicate or sensitive skin.

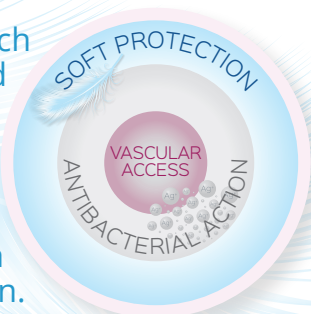


REFERENCES

Data and technical specifications can be changed without notice.

CODE	PLASTER (cm)	PAD (cm)	HOLE ø (cm)	PCS / PACK	PACKS / BOX
S0810/AG	10 x 8	5 x 4	0.4	1	50

Traditional vs. SoftTouch adhesives


Traditional adhesives

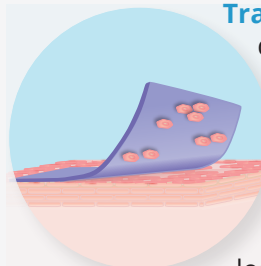

do not adhere evenly to the skin; instead, they leave interspaces on the surface which are subsequently filled. A traction occurs on the skin that makes the fix stronger and more suitable for long-term treatment.

SoftTouch adhesives

adhere evenly to the skin and have low surface tension. Removing them is less traumatic, making this kind of adhesive ideal for frequent treatment or for delicate and sensitive skin.

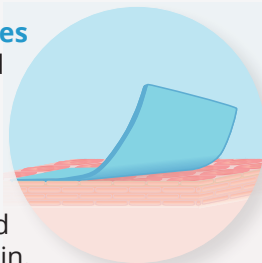

INSTRUCTION FOR USE

- Disinfect the exit-site avoiding iodine-based disinfectants and let it dry completely.
- Lift the two flaps of the protective film of the plaster.
- Place the pad in contact with the skin at the point of exit of the catheter and fix it by gentle pressure on the surrounding plaster.
- Complete the removal of the protective film.
- Maintain **EXIT-PRO AG SoftTouch** on site up to one week, if the exit-site does not show initial sign of infection.



WARNINGS

- Do not use in patients with hypersensitivity to silver.
- Do not use during examinations such as X-rays, ultrasound, diathermy, magnetic resonance or radiotherapy treatments.
- Do not use in combination with iodine-based solutions.
- Apply on dry skin.
- It does not replace the appropriate systemic therapies necessary in presence of clinical infections.
- Read the instructions for use carefully before application.

CLASSIFICATION • CLASS III (RULE 13 MDD 93/42 EEC)

BIOCOMPATIBLE • ACCORDING TO UNI EN ISO 10993 STANDARDS

STERILIZATION • ETHYLENE OXIDE

EXPIRATION • 3 YEARS



• SINGLE-USE